SVF Flow Controls

Series SB7 Ball Valve

Three-Piece CleanFLOW[™] High Purity Ball Valve Sizes 1/2" ~ 6"

CleanFLOW[™] SB7 ball valves are engineered to be a true process piping component to specifically meet the demanding processes found in the pharmaceutical and food & beverage industries. The "Tube-ID" port opening is dimensionally identical to the adjacent tubing to comply with ASME-BPE guidelines. All materials comply with FDA, USDA and 3A requirements.

✓ ASME-BPE compliant ✓ ASME-BPE compliant ✓ Silicone Free ✓ Cavity filled TFM1600[™] seat option available ✓ Complete 316L Stainless Steel cast construction ✓ Drainable design with "Tube-ID" dimensions ✓ Weld bosses for easy purge porting on ends

- \checkmark ISO 5211 mounting pad for easy actuation
- Encapsulated body seals to facilitate welding without disassembly
- End connections include Tri-Clamp and Extended Tube O.D.
- ✓ Controlled delta ferrite chemistry
- Chevron (V-Ring) Stem Seals assures low friction and leak tight performance (Sizes 2-1/2" ~ 4" only)
- \checkmark Standard interior finish is 15Ra or better
- 🖌 ETO ends are designed for Orbital Welding
- \checkmark Exclusive "Fine Adjust" handle for precise positioning on sizes 1/2" ~ 2"
- \checkmark Locking device is standard on 3" and 4" valves


MATERIALS OF CONSTRUCTION


	J.						
ITEM #	DESCRIPTION	MATERIALS SPECIFICATIONS (Additional options available)					
1	Body	316L Stainless Steel (ASTM A351 CF3MN					
2	End Connector	316L Stainless Steel (ASTM A351 CF3MN)					
3	Ball	316L Stainless Steel (ASTM A351 CF3MN)					
4	Stem	316L Stainless Steel (ASTM A351 CF3MN)					
5	Seat	TFM1600™					
6	Body Seal	PTFE/TFM1600™					
7	Stem Seal	em Seal TFM1600™					
8	Locking Device (Optional)	304 Stainless Steel					

What do you need today?™





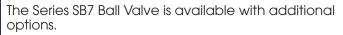

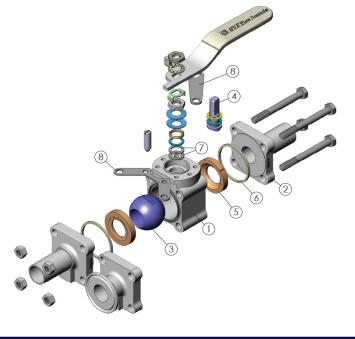

There are no elastomers used for any of the components on the SB7 valves. Only PTFE and TFM are used for seats and seals.





SVF Flow Controls • 5595 Fresca Drive • La Palma, CA 90623 • Tel: 800.783.7836 • FAX: 562.802.3114 • Sales@SVF.net Specifications subject to change w/o notice. All Data Sheets on our website supersede prior publications • (SVF.SB7.1118)

SVF Flow Controls

Series SB7 Ball Valve


Three-Piece CleanFLOW[™] High Purity Ball Valve


Sizes 1/2" ~ 6"

DIMENSIONS, WEIGHT, CV, TORQUE

Size	A-1		A	A-2		В		D		L		W		Weight		Toro	ue*
	in.	mm	in.	mm	in.	mm	in.	mm	in.	mm	in.	mm	lbs	kg		in-Ibf	Nm
1/2″	5.50	140	3.50	89	0.37	9	5	114	2.22	56	0.065	1.7	2	0.9	8	60	7
3/4″	6.00	152	4.00	102	0.62	16	5	114	2.28	58	0.065	1.7	2	0.9	29	60	7
" ا	6.50	165	4.50	114	0.87	22	6	146	2.56	65	0.065	1.7	4	1.8	66	100	11
1-1/2″	7.50	191	5.50	140	1.37	35	7	178	3.09	78	0.065	1.7	8	3.6	192	200	23
2″	8.00	203	6.25	159	1.87	47	7	178	3.45	88	0.065	1.7	13	5.9	434	250	28
2-1/2″	9.50	241	6.75	171	2.37	60	10	254	5.35	136	0.065	1.7	23	10.4	723	450	51
3″	10.50	267	7.00	178	2.87	73	14	348	6.55	166	0.065	1.7	31	14.1	1124	1300	147
4″	12.50	318	8.50	216	3.83	97	22	559	7.14	181	0.083	2.1	46	20.9	2100	1400	158
6″	16.00	406	17.00	432	5.78	147	26	660	12.00	305	0.109	2.8	196	88.5	4700	4160	470

SB7 - PRESSURE/TEMPERATURE CHART

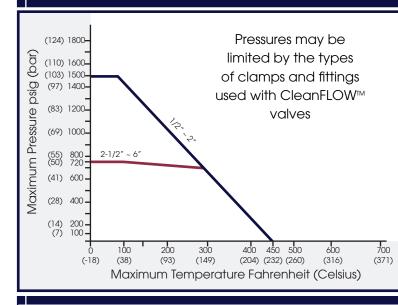

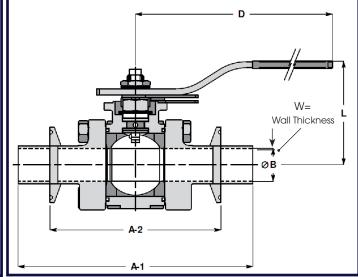

* At full differential pressure for clean fluids with TFM1600™

Non-Cavity Filled Seats

HOW TO ORDER SERIES SB7 BALL VALVES

Please refer to our website for our comprehensive How to Order Guide for Series SB7 Ball Valves. SVF Flow Controls

Series SB7 Ball Valve

Three-Piece CleanFLOW[™] High Purity Ball Valve

Sizes 1/2" ~ 6"

Ordering Code Sequence (Columns 1 thru 11)											
1	2		3	4	5		6				
			ENDS hless Steel 51 CF3MN	BALL 6 = 316L Stainless Steel ASTM A351 CF3M	STEM 6 = 316L Stainless Steel ASTM A351 CF3M		SEAT MATERIAL A = TFM1600™ Q = TFM1600™ Cavity Filled				
7	8		9	10			11				
BODY SEAL T = PTFE A = TFM1600™ (6" Size Only)	BODY SEALEND CONNECTIO= PTFETR0 = Tri-Clamp Ends= TFM1600™ETO =		VALVE SIZE 005 = 1/2" 007 = 3/4" 010 = 1" 015 = 1-1/2" 020 = 2" 025 = 2-1/2" 030 = 3" 040 = 4" 060 = 6"	OPTIONS* 00 = Standard Lever LK = Locking Device AA = Oval Handle & Locking Device AC = Locking Device Anti-Static Device AG = ISO Stem Extensic & Locking Device JA = Oval Handle, ISO Stem Extensic & Anti-Static Dev JE = Locking Device ISO Stem Extensic & Anti-Static Dev JE = Locking Device	e & ce sion re on vice , on	00 = None XC = Oxyg EP = Electro SB = 10Ra I AB = Electro 10Ra ID AD = Oxyg Electrop JB = Electro	D Finish opolished & Finish en Cleaned & polished polished, Cleaned &				

